



# **Practical Course Of Physical chemistry**

**Prepared by**

**M.S. Kosrat and Karzan**

### **Contents:**

**Experiment (1):** Variation of conductance with concentration (strong & weak) electrolytes.

**Experiment (2) :** Determine Hydrolysis constant by a conductimetric method.

**Experiment (3):** Hydrolysis constant by a pH method.

**Experiment (4):** Dissociation constant of the weak acid (Conductometric method).

**Experiment (5):** Conductometric titration of strong, weak mixture acid with strong base.

**Experiment (6):** Solubility and solubility product from conductivity measurements.

**Experiment (7):** Dissociation constant by the pH method.

### Some basic concepts

In experiment of chemistry, mostly students require standard solution to study physical properties.

#### Standard solution:

A solution of known concentration is called standard solution. It contains a known weight of substance in a definite volume of solution. Generally concentrations are expressed in terms of

- i) **Normality:** It is the number of gram equivalents of solute present in one liter of solution.

$$\text{Normality (N)} = \frac{\text{Number of gram equivalents}}{\text{Liter of solution}}$$

- ii) **Molarity:** It is the number of moles of solute present in one liter of solution.

$$\text{Molarity (M)} = \frac{\text{Number of moles}}{\text{Liter of solution}}$$

- iii) **Molality:** It is the number of moles of solute dissolved in 1000 gms of solvent.

$$\text{Molality (m)} = \frac{\text{Number of moles of solute}}{\text{Wt. of solvent in kg (1000 gms)}}$$

- iv) **Mole fraction:** The mole fraction of a component in a solution is defined as the number of moles of that component divided by the total number of moles of all components in the solution. If a solute is dissolved in solvent then

Moles of solute

The mole fraction of solute = \_\_\_\_\_

Moles of solute + Moles of solvent

Moles of solvent

The mole fraction of solvent = \_\_\_\_\_

Moles of solute + Moles of solvent

- v) **Volume fraction:** If two liquids are miscible to form a solution then the volume fraction can be determined as follows:

 $v_1$ Volume fraction of one ( $V_1$ ) = \_\_\_\_\_ $v_1 + v_2$  $v_2$ Volume fraction of other ( $V_2$ ) = \_\_\_\_\_ $v_1 + v_2$ 

- vi) **Weight fraction:** It is the ratio of weight of solute to the total weight of solution.

Wt. of solute

Wt. fraction of solute = \_\_\_\_\_

Total wt. of solution

Wt. of solvent

Wt. fraction of solvent = \_\_\_\_\_

Total wt. of solution

If weight fraction multiplied by 100 it is termed as percentage by weight.

vii) **Percentage weight of solute by volume of solvent:** It is the weight of solute dissolved in 100 ml of solution: (Wt/V) and (Wt/Wt). also there is percentage volume of solute by volume of solvent (V/V).

viii) **Parts per thousand (P.P.T.):** This is a method of expression concentration of solute if solution is sufficiently dilute.

If 20 gm of can sugar is dissolved in 1000 gm of water (or 1000 ml) then the resultant solution is said to be 20P.P.T. it can be expressed as

$$\text{P.P.T} = \text{gm}/1000 \text{ gm of solution.}$$

ix) **Parts per thousand (P.P.m.):** This is the best method of express concentration of solute when solution is expected to be very dilute.

If 20 gm of can sugar is dissolved in one million gm of water (or 1000 ml) then the resultant solution is called 20 P.P.m. it can be expressed as

$$\text{P.P.m} = \text{gm}/1000000 \text{ gm of solution.}$$

x) **Strength of solution:** The strength of solution is always expressed in terms of grams per liter.

If the concentration is given in terms of normality then

$$\text{Strength} = \text{Normality} \times \text{Equivalent wt.}$$

And if the concentration is given in terms of normality then

$$\text{Strength} = \text{Molarity} \times \text{Molecular wt.}$$

xi) **Dilution:** A solution is diluted when extra solvent is added. The concentration of dilute solution can be calculated by using this equation which known as Dilution law

$$M_1 V_1 = M_2 V_2 \text{ Or } N_1 V_1 = N_2 V_2.$$

**Preparation of solution****1) Preparation solution from solid material**

$$\text{Wt}_{(\text{gm})} = \frac{M \times V_{\text{ml}} \times \text{M.wt}}{1000}$$

**For example:**

prepare 0.1 M of NaOH in 100 ml of distilled water.

M.wt of NaOH = 23+16+1= 40 gm/mol

V = 100 ml

M = 0.1 M

Wt = ?

$$\text{Wt}_{(\text{gm})} = \frac{M \times V_{\text{ml}} \times \text{M.wt}}{1000} = \frac{0.1 \times 100 \times 40}{1000} = 0.4 \text{ gm}$$

Take 0.4 gm of NaOH and dissolved in quantity of water, when it dissolved completely add it into 100 ml standard flask and fill it up to the mark.

**For example:**

Prepare 0.1 N of  $\text{Ca(OH)}_2$  in 100 ml of distilled water.

Eq.wt of  $\text{Ca(OH)}_2 = \text{M.wt} / 2 = 74/2 = 37 \text{ gm / mol}$  , Eq.wt = M.wt /  $n_{(\text{no. of H}^+, \text{OH}^- \text{ and e}^-)}$

$V = 100 \text{ ml}$

$N = 0.1 \text{ N}$

$\text{Wt} = ?$

$$\text{Wt}_{(\text{gm})} = \frac{N \times V_{\text{ml}} \times \text{Eq.wt}}{1000} = \frac{0.1 \times 100 \times 37}{1000} = 0.37 \text{ gm}$$

Take 0.37 gm of  $\text{Ca(OH)}_2$  and dissolved in quantity of water, when it dissolved completely add it into 100 ml standard flask and fill it up to the mark.

**2) Preparation solution from liquid material**

For preparation of dilute solution from concentrate liquid can be use these information. Density, % and M.wt or Eq.wt.

$$\text{Sp.gr} = \frac{\text{Density of liquid}}{\text{Density of water at } 4^\circ\text{C}}$$

Density of water at  $4^\circ\text{C} = 1$ , i.e.  $\text{Sp.gr} = \text{Density}$

% is the percentage of concentrate liquid in bottle.

$$M = \frac{\text{Sp.gr} \times \% \times 1000}{\text{M.wt}} \quad \text{Or} \quad N = \frac{\text{Sp.gr} \times \% \times 1000}{\text{Eq.wt}}$$

By using dilution law:

$$M_1 V_1 = M_2 V_2 \quad \text{Or} \quad N_1 V_1 = N_2 V_2$$

**For example:**

Prepare 0.5 M and N of  $\text{H}_2\text{SO}_4$  in 100 ml of distilled water. If you know density of concentrate  $\text{H}_2\text{SO}_4 = 1.84 \text{ gm/ml}$ , M.wt= 98 gm/mol and % = 96.

$$M = \frac{\text{Sp.gr} \times \% \times 1000}{\text{M.wt}} \quad \text{for} \quad N = \frac{\text{Sp.gr} \times \% \times 1000}{\text{Eq.wt}}$$

$$= \frac{1.84 \times 96/100 \times 1000}{98} \quad \text{for} \quad = \frac{1.84 \times 96/100 \times 1000}{98/2}$$



$$= \frac{1.84 \times 0.96 \times 1000}{98} \text{ for } = \frac{1.84 \times 0.96 \times 1000}{49}$$

$$M = 18.02 \quad \text{for} \quad N = 36.04$$

Now, for Prepare 0.5 M and N of  $\text{H}_2\text{SO}_4$  in 100 ml of distilled water, using dilution law.

$$M_1 V_{1(\text{conc})} = M_2 V_{2(\text{dilute})} \quad \text{for} \quad N_1 V_{1(\text{conc})} = N_2 V_{2(\text{dilute})}$$

$$18.02 \times V_{1(\text{conc})} = 0.5 \times 100 \quad \text{for} \quad 36.04 \times V_{1(\text{conc})} = 0.5 \times 100$$

$$V_{1(\text{conc})} = 2.77 \text{ ml} \quad V_{1(\text{conc})} = 1.38 \text{ ml}$$

Take calculated volume of concentrate liquid and put it in standard flask fill up to the mark by distilled water.

## Electrochemistry

### Strong electrolytes:

According to Arrhenius theory, strong electrolytes are completely dissociated at any concentration. It means if a strong electrolyte is dissolved in water, it will ionize completely into cations and anions. In other words it can be stated that the ionization is 100% and question of calculating degree of dissociation and dissociation constant, does not arise at all. Therefore, Ostwald's law of dilution cannot be made applicable to the strong electrolytes.

If concentrated solution of strong electrolyte is diluted the ions get more space to movement and therefore conductivity of strong electrolyte increases with dilution up to certain extent and attains limiting value.

Electrolytes like HCl, NaOH, KCl have high value of conductance even at low dilution.

### Weak electrolytes:

Weak electrolytes are partially dissociated in concentrated solution. The extent of dissociation is called degree of dissociation. The degree of dissociation increases with dilution and therefore it can be stated that conductance of weak electrolyte increases with dilution. The degree of dissociation of weak electrolyte needs to be calculated at particular concentration and hence dissociation constant can be calculated by applying Ostwald's law of dilution.

If degree dissociation of an electrolyte is considered to be 0.2, it means 20% of electrolyte is dissociated into ions and 80% of electrolyte will be in undissociated condition. Therefore, only 20% of electrolyte will give conductance value. The remaining 80% of electrolyte will not take part in conducting the electricity. As the conductivity value is proportional to number of dissociated molecules the degree of dissociation can be calculated by conductivity measurements.

Weak electrolytes like acetic acid, ammonium hydroxide, sodium acetate and ammonium chloride have low value of conductance at low dilution as they are partially ionized in the solution.

**Observed conductivity:**

The conductivity recorded by any conductivity cell is termed as observed conductivity.

**Specific conductivity:**

The specific conductivity of an electrolytic solution may be define as the conductance of a solution is placed between electrodes 1 cm apart and having cross-sectional area 1 sq.cm of each electrode. Specific conductivity means conductance of one cubic centimeter of electrolyte. The value of specific conductance depends upon the number of ions (concentration of ions) present in one cubic centimetre of solution.

The specific conductivity decreases with increasing dilution because number of ions per unite volume decreases.

**The specific conductivity = cell constant x observed conductance**

**Equivalent conductivity (  $\kappa_e$  ):**

It is the power of all ions produced by one gram equivalent of an electrolyte at any particular concentration.

$$\kappa_e = K_c \frac{1000}{C}$$

Where,  $\kappa_e$  is equivalent conductivity,  $K_c$  is specific conductivity and C is the concentration.

The equivalent conductivity increases with increasing dilution and attains a limiting value at infinite dilution. The limiting value is called equivalent conductance at infinite dilution, denoted

by . The increase of equivalent conductivity is due to the increase of ions as a result of dissociation. As the ions get sufficient space to move, migration velocity increases.

**Molecular conductivity:**

Molecular conductance of a solution is the conductivity of all the ions produced by dissociation of one gram molecular weight of electrolyte when dissolved in Vml.

$$\mu_v = K_c \times 1000 / M$$

where,  $\mu_v$  is the molecular conductivity.

**Ionic conductance:**

Weak electrolytes are partially ionized in solution and ionization increases with dilution. The conductivity is due to mobility of cations and anions. Hydrogen ions ( $H^+$ ) is having high mobility than  $Na^+$  and  $CHCOO^-$  ions because it is lighter than sodium ion and acetate ion. It is not possible to calculate of weak electrolyte experimentally. Then can be calculated by applying Kohlraush law.

$$= \text{Ionic conductivity of cation} + \text{Ionic conductivity of anion}$$

$$= \lambda_c + \lambda_a$$

The molar ionic conductivity values are referred from literature.

Molar/Ionic conductivity values as 25 °C

Cation	$\lambda^\circ$ (S-cm <sup>2</sup> /mol)	Anion	$\lambda^\circ$ (S-cm <sup>2</sup> /mol)
H <sup>+</sup>	349.6	OH <sup>-</sup>	199.1
Li <sup>+</sup>	38.7	F <sup>-</sup>	55.4
Na <sup>+</sup>	50.10	Cl <sup>-</sup>	76.35
K <sup>+</sup>	73.50	Br <sup>-</sup>	78.1
Rb <sup>+</sup>	77.8	I <sup>-</sup>	76.8
Cs <sup>+</sup>	77.2	NO <sub>2</sub> <sup>-</sup>	71.8
Ag <sup>+</sup>	61.9	NO <sub>3</sub> <sup>-</sup>	71.46
NH <sub>4</sub> <sup>+</sup>	73.5	ClO <sub>3</sub> <sup>-</sup>	64.6
Ethylammonium	47.2	ClO <sub>4</sub> <sup>-</sup>	67.3
Diethylammonium	42.0	IO <sub>4</sub> <sup>-</sup>	54.5
Triethylammonium	34.3	HCO <sub>3</sub> <sup>-</sup>	44.5
Tetraethylammonium	32.6	H <sub>2</sub> PO <sub>4</sub> <sup>-</sup>	57
Tetra- <i>n</i> -butylammonium	19.5	HSO <sub>3</sub> <sup>-</sup>	50
Dimethylammonium	51.8	HSO <sub>4</sub> <sup>-</sup>	50
Trimethylammonium	47.2	HC <sub>2</sub> O <sub>4</sub> <sup>-</sup>	40.2
Tetramethylammonium	44.9	HCOO <sup>-</sup>	54.6
Piperidinium	37.2	CH <sub>3</sub> COO <sup>-</sup>	40.9
		C <sub>6</sub> H <sub>5</sub> COO <sup>-</sup>	32.4

**Exp. 1****Variation of conductance with concentration (strong & weak ) electrolytes****Theory**

The specific conductance  $\cdot K$  of any conductor is defined as the (reciprocal of the specific resistance) that is

$$K = (L/aR) \quad (\text{ohm}^{-1} \text{cm}^{-1})$$

Where  $R$  is the resistance,  $L$  is the length in centimeters and  $a$  is the cross sectional area in square centimeters of the conductor. The equivalent conductance by multiplying the specific conductance by the volume of the solution measured in cubic centimeters containing 1 g equivalent of electrolyte

$$= K \cdot v \quad (\text{ohm}^{-1} \text{cm}^2 \text{eqt.}^{-1})$$

Hence, if  $c$  is the concentration of the solution in gram equivalents per liter,  $v$  is equal to  $1000/c$ , and

$$= K \cdot 1000/c \quad (\text{ohm}^{-1} \text{cm}^2 \text{eqt.}^{-1})$$

Note that an equivalent is the associated with the transfer of one faraday of electricity. Hence, if one mole carries  $n$  faradays of electricity there will be  $n$  equivalent in one mole, Only in the case where one mole is associated with the transfer of one faraday the mole and equivalent are identical.

Both specific and equivalent conductance of a solution vary with concentration .Specific conductance increases with increase in concentration .While equivalent conductance increases with dilution, a limiting value  $\kappa_0$  being obtained at infinite dilution (zero concentration). The equivalent conductance at zero concentration of an electrolyte is not the same as the equivalent conductance of the pure solvent .it is the equivalent conductance of the solution when the ions of the electrolyte

are on infinite distance apart, and hence free from interaction.

For certain electrolytes, Kohlrausch suggested the empirical relationship

$$\kappa = \kappa_0 - b \cdot c$$

Where  $b$  is a constant .These electrolytes are known as strong electrolytes and a plot of  $\kappa$  vs.  $c$  is linear for concentration up to 0.01 N .The equation is accurate to within a few percent ,however ,for concentrations of up to 0.1 N in the case of unit-univalent electrolytes .dilute aqueous solutions Of many other electrolytes ,e.g. acetic acid ,do not give a linear graph .For

electrolytes, known as weak electrolytes, the value of the equivalent conductance increases rapidly as  $c$  approaches zero, and the graph cannot be accurately extrapolated to obtain a value of  $\Lambda^0$ .



### Procedure

1. Prepare 50 ml of 0.1 M (M/10) standard potassium chloride solution to be used for determining cell constant. The electrodes are washed with distilled water and the conductance of it 50ml recorded. And read conductance of KCl for eight diluted (0.01, 0.02, 0.03, 0.04, 0.05, 0.06, 0.07 and 0.08 molar solutions which prepared from 0.1M potassium chloride stock solution).
2. Prepare 50 ml of 0.1 M acetic acid. Read the conductance of 0.1 M acetic acid and for eight diluted (0.01, 0.02, 0.03, 0.04, 0.05, 0.06, 0.07 and 0.08 molar solutions which prepared from 0.1M acetic acid stock solution).

### Calculation

1. Calculate the cell constant.  $\kappa = G K \text{ cm}^{-1}$ , where **K** is the “cell constant” and G is the observed conductance, k is the specific conductivity for M/10 KCl at 25°C = 0.01288 mhos.
2. Equivalent conductance (  $\Lambda$  ),  


$$= 149.82 - 93.85 \sqrt{c} + 94.9c(1 - 0.2274 \sqrt{c})$$
3. Specific conductance (k),  $K = \frac{\kappa}{C} \text{ (cm}^{-1}\text{)}$
4. Molar conductance for strong and weak Electrolyte,  $m = k/C$ , where k is the specific constant.
5. Plote  $\Lambda$  vs.  $\sqrt{c}$  thane from the obtained curve show the type of week electrolyte.

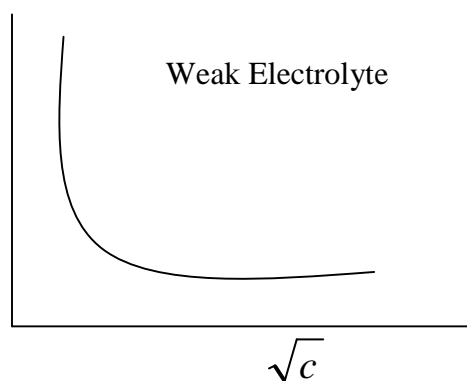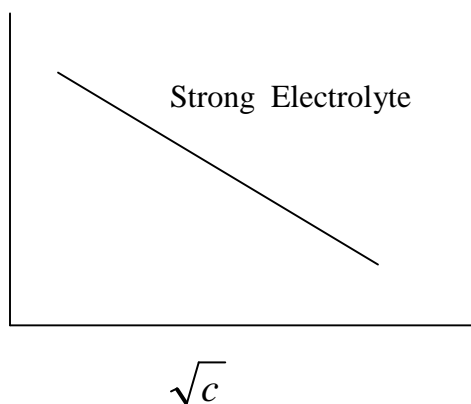



**Exp.2****Hydrolysis constant by a conductimetric method****Theory**

When the salt of a strong acid and weak base such as aniline chloride is dissolved in water, hydrolysis occurs, the cations of the salt interacting with the solvent to form the conjugate base of the salt and the acid of the solvent.

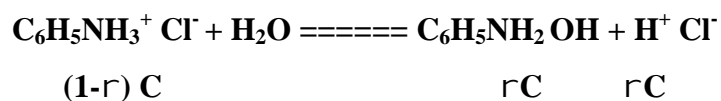

If fraction  $\alpha$  of the salt hydrolysis, the hydrolysis constant is given by :-

$$K_h = r^2 C / (1-r)$$

Where  $v$  is the volume in liters containing one mole of salt. the conductivity of the salt in which hydrolysis has occurred will be made up partly of the conductivity due to the ions produced by hydrolysis.

If excess of aniline is added in this solution the hydrolysis solution is suppressed and therefore conductivity mainly depends upon salt solution in presence of aniline and strong acid HCl formed during the hydrolysis.

For determination of degree of hydrolysis of the salt by conductivity method, following different measurements of conductivity are necessary.

- i) The conductivity of aniline hydrochloride in water at certain concentration.
- ii) The conductance of salt solution in presence of base i.e aniline.
- iii) The conductivity of strong acid formed during the hydrolysis process.

Let  $r$  be the degree of hydrolysis. Therefore,

$$\kappa = (1-r) \kappa' + r \kappa''$$

$$\kappa - \kappa' = r (\kappa'' - \kappa')$$

$$r = (\kappa - \kappa') / (\kappa'' - \kappa')$$

Where  $\kappa$  is the equivalent conductivity of aniline hydrochloride salt in water,  $\kappa'$  is the equivalent conductivity of aniline hydrochloride in presence of aniline and  $\kappa''$  is the equivalent conductivity of HCl at that concentration.

Hydrolysis constant can be calculated by Ostwald's law of dilution.

$$K_h = r^2 C / (1-r)$$

And it is possible to calculate the affinity constant of aniline  $K_b$  from the relationship .

$$K_h = (K_w / K_b)$$

### Procedure

1. Prepare 50 ml of 0.1 N (N/10) standard potassium chloride solution to be used for determining cell constant. The electrodes are washed with distilled water and the conductance of it 50ml recorded.
2. Prepare N/32 solution of aniline chloride in water take 50ml of this solution is placed in the conductivity cell and note down the conductivity at 25°C.
3. Prepare N/64 and N/128 solution from N/32 solution of aniline chloride by dilution, and read their conductivity.
4. Prepare also N/32, N/64 and N/128 solution of aniline chloride by using N/32 solution of aniline as a solvent and read their conductivity at 25°C.

### Calculation

1. Calculate the cell constant. Specific constant = cell constant of KCl x observed conductance, the specific conductivity for M/10 KCl at 25°C = 0.01288 mhos.
2. Specific constant ( $K_c$ ) for each solution.
3. Equivalent conductance for each solution by this relation,  $\kappa = K_c \cdot 1000/C$ .
4. Degree of hydrolysis at each concentration,  $r = \frac{\kappa - \kappa_0}{\kappa - \kappa_0'}$ .
5. Hydrolysis constant of aniline chloride,  $K_h = r^2 C / (1-r)$ .
6. Dissociation constant of aniline,  $K_a = (K_w / K_b)$ .

$\kappa$  for hydrochloric acid at 25°C

N/32    393     $\text{cm}^2 \text{eqt.}^{-1}$

N/64    399     $\text{cm}^2 \text{eqt.}^{-1}$


N/128    401     $\text{cm}^2 \text{eqt.}^{-1}$

$K_w$  for water at 25 °C =  $1.01 \times 10^{-14} (\text{g ions})^2 \text{L}^{-2}$

**Exp.3****Hydrolysis constant by a pH method****Theory**

The value of hydrogen concentrated predicate from hydrolysis of strong acid and weak base or from strong base and weak acid depend of the formation of pH or from acid value of solution:

$$\text{pH} = -\log[\text{H}^+]$$

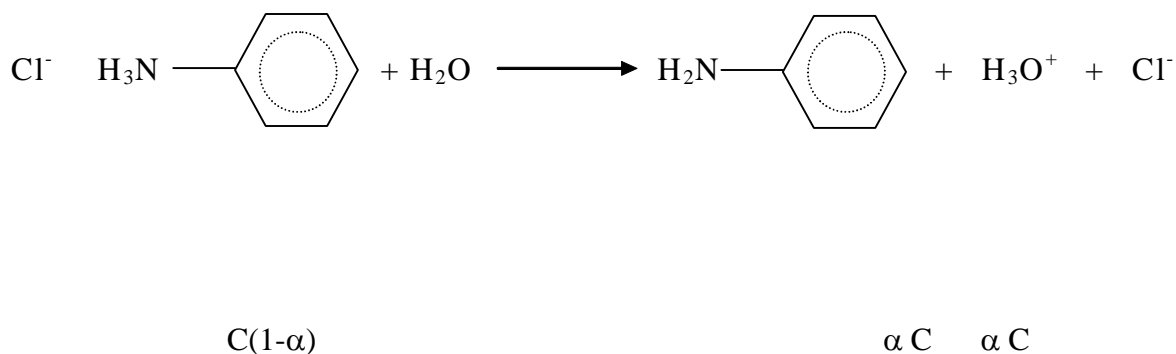

Where  $\alpha$  is degree of hydrolysis, and C concentration of solution.

$$K_h = (r^2 C) / (1-r)$$

$K_h$  is hydrolysis constant

**Procedure**

1. Prepare 50 ml of 0.1M of aniline hydrochloride and dilute from it 0.001 , 0.0001 and 0.00001 M .
2. Read pH of solutions after 20 min, when at hydrolysis is complete.

**Calculation**

1. calculation the hydrogen ion concentration at each stoichiometric concentration by

$$\text{pH} = -\log[\text{H}^+]$$

2. calculate degree of hydrolysis of each concentrated ion by

$$[\text{H}^+] = r \cdot C$$

3. calculate the hydrolysis constant of each concentration by


$$K_h = (r^2 C) / (1-r)$$

4. calculate ionic hydrolysis of aniline by

$$K_b = K_w / K_h$$

**Exp. 4****Dissociation constant of the weak acid (Conductometric method)****Theory**

According to the so-called Ostwald dilution law, for a binary weak electrolyte, e.g. the weak acid HA.

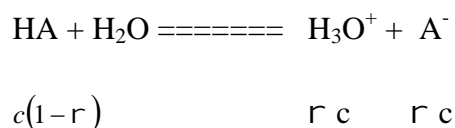

$$K_d = \frac{\alpha^2 c}{(1-\alpha)}$$

Where  $\alpha$  is the degree of dissociation,  $c$  is the concentration, and  $K_d$  is a constant sometimes referred to as the classical dissociation constant according to the Arrhenius dissociation theory, the degree of dissociation is given by the ratio of the equivalent conductance to the conductance at infinite dilution.

$$\alpha = \Lambda / \Lambda_o$$

Hence, substituting for  $\alpha$  is the expression for the dissociation constant we have.

$$K_d = \frac{\alpha^2 \cdot c}{\Lambda_o(\Lambda_o - \Lambda)}$$

$$\Lambda^2 c = K_d \Lambda_o^2 - K_d \Lambda_o \Lambda$$

$$\Lambda \cdot c = K_d \Lambda_o^2 \cdot \frac{1}{\Lambda} - K_d \Lambda_o$$

$$\Lambda \cdot c = K_d \left( \frac{\Lambda_o^2 - \Lambda_o}{\Lambda} \right)$$

A plot of  $\Lambda \cdot c$  vs.  $1/\Lambda$  will thus give a straight line as shown in Fig.1 at

$C = 0$ ,  $OX = 1/\Lambda_o$ , that is  $\Lambda_o = 1/OX$ . at  $1/\Lambda = 0$ ,  $OY = -K_d \cdot \Lambda_o$ , therefore  $(OX)(OY) = -K_d$


Where  $K_d$  is the classical dissociation constant

### Procedure

1. Prepare 50 ml of 0.1 M (M/10) standard potassium chloride solution to be used for determining cell constant. The electrodes are washed with distilled water and the conductance of it 50ml recorded.
2. Prepare (0.01, 0.02, 0.03, 0.04, 0.05, 0.06, 0.08) solution from 0.1M of acetic acid stock solution. , and read their conductivity.

**Calculation**

1. Calculate the cell constant. Specific constant = cell constant of KCl  $\times$  observed conductance, the specific conductivity for M/10 KCl at 25°C = 0.01288 mhos.
2. Specific constant ( $K_c$ ) for each solution.
3. Equivalent conductance for each solution by this relation,  $\kappa_c = K_c \cdot 1000/C$ .
4. Plot of  $\Lambda_c$  vs.  $1/\Lambda$ .

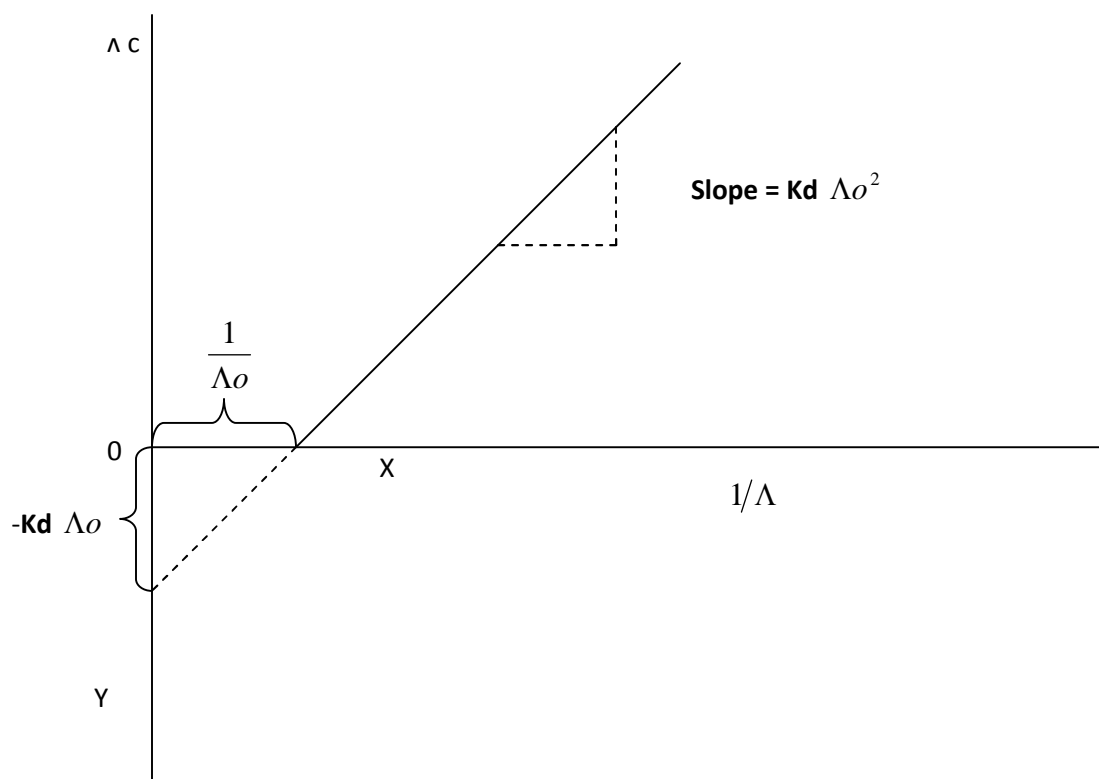



**Exp.5****Conductometric titration of strong, weak mixture acid with strong base****Theory**

On the addition of a solution of a strong acid, there is initially a decrease in conductivity due to the replacement of the fast moving hydrogen ion by the slower alkali metal ion the conductivity subsequently increases again with the increase in concentration of hydroxyl ions when the titration is complete. Hence if the conductivity is plotted against the volume of base solution added the end point of the titration given by a Sharp inflection in the graph.

If strong base, however, is added to a solution of a weak acid, there is an initial gradual increase in the conductivity as the feebly dissociated weak electrolyte is converted into its salt, a stronger electrolyte. Subsequently at the end point a much more rapid conductivity increase occurs due to the presence of excess hydroxyl ions. If a mixture of strong acid and a weak acid is titrated by a strong base, an initial decrease in conductivity is followed by an increase, and then by a more rapid increase thereby giving two end points, see Figs.




**Procedure**

1. Prepare 0.1 M from each sodium hydroxide, hydrochloric acid and acetic acid in 50 ml.
2. Dilute 10 ml from hydrochloric acid in 50 ml and immersed in beaker to measured conductivity the 0.1 M sodium hydroxide solution is then run to acid solution and the conductance of the solution is measured after addition of 2,4,6,8,10 ..... ml of sodium hydroxide to after end point is present before the measurement of the conductance is made the solution should be thoroughly stirred.
3. Repeat this experiment using acetic acid instead of hydrochloric acid.
4. Prepare mixture of 5ml hydrochloric acid and 5ml acetic acid and diluted to 50 ml distilled water and repeat adding of sodium hydroxide and measure conductivity

### Calculation

1. Plotted graph conductivity vs. volume of sodium hydroxide for each acid and mixture.
2. Obtain from plotted the end point for each acid and mixture.
3. Determine the concentration of each acid using

$$(N_1 \cdot V_1)_{\text{of NaOH}} = (N_2 \cdot V_2)_{\text{of acid or mixture}}$$

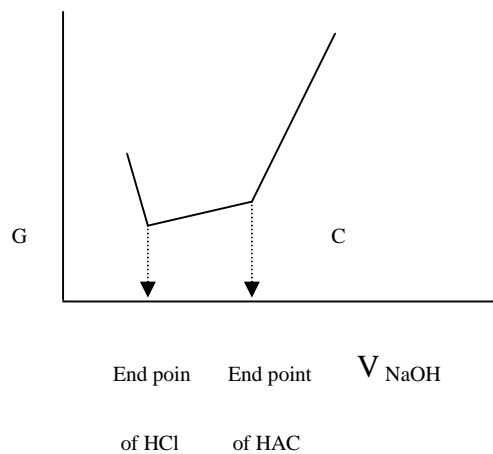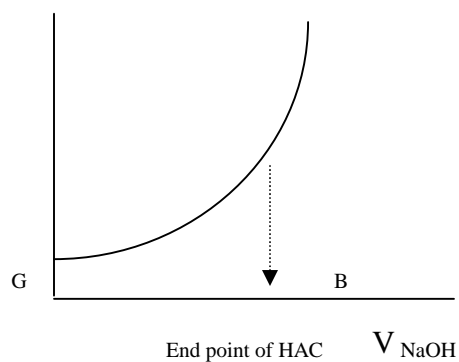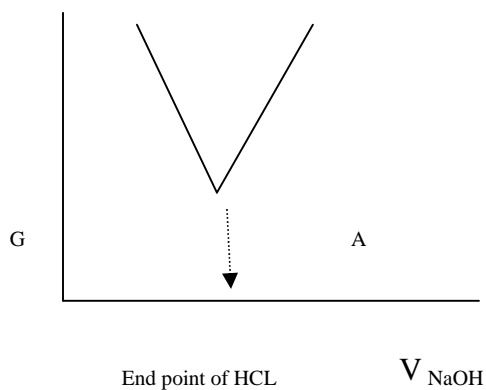

**Exp.6****Solubility and solubility product from conductivity measurements****Theory**

The Solubility product of an electrolyte may be determined from conductivity measurements. The conductance method is application for finding the concentration of the saturated solution provided that the salt is not hydrolysis and that solubility is not too high. It may then be assumed that the ions possess their limiting conductivities these must of course be known for the same temperature, measurement of the specific conductivity of the saturated solution leads to a value for the concentration.

**Procedure**

1. Prepare 50 ml of 0.1 M (M/10) standard potassium chloride solution to be used for determining cell constant. The electrodes are washed with distilled water and the conductance of it 50ml recorded.
2. Prepare silver chloride, silver bromide and silver iodide by preparation 0.1M from each potassium chloride, potassium bromide and potassium iodide in 50ml and prepare 0.1M of silver nitrate in 250 ml after this prepare mixed each 50ml of three sparingly soluble salts with 50ml from 0.1M silver nitrate, after this filtered each solution to measure conductivity.

**Calculation**

1. Calculate the cell constant. Specific constant = cell constant of KCl x observed conductance of 0.1 M KCl, the specific conductivity for M/10 KCl at 25°C = 0.01288 mhos.
2. Specific constant ( $K_c$ ) for each solution.
3. Equivalent conductance for each solution by this relation,  $\kappa_c = K_c \cdot 1000/C$ .
4. Calculate equivalent conductance for saturated solution by,  $\kappa_c = \kappa + \kappa_a$ .
5. Calculate solubility by,  $s = \kappa / \kappa_c \times 1000$ .
6. Calculate the solubility product by,  $K_{sp} = S^2$

Not /  $\kappa_c$  and  $\kappa_a$  it is equivalent ionic conductivities at finite dilution can be determined from references.

**Exp. 7****Dissociation constant by the pH method****(pH titration curve)****Theory**

One form of the glass electrode consists of a glass tube terminating in a thin walled glass bulb containing a platinum contact immersed in a buffer solution. Its potential (EG) when immersed in a buffer solution is dependent on the pH of the same test solution and is given by the equation

$$EG = \text{constant} + \frac{RT}{F} \ln aH^+$$

Where  $aH^+$  is the activity of the hydrogen ions in the test solution. The precise potential of this electrode is not of great interest as it includes an asymmetry potential of this electrode is not of great interest as it includes an asymmetry potential (possibly due to strain at the glass solution interface)

This exists across the membrane even when the two solutions on either side of the membrane have identical activity values. In practice, however assembly shown below. Consequently this pH can be used as an indicator in acid – alkali titrations, there being a relatively large change in pH at the end point.

In this experiment the glass electrode is connected in series with saturated dip-type calomel electrode as shown:

Solution of constant pH | glass membrane | solution of unknown pH | Calomel electrode

The change in pH during the titration of sodium hydroxide with weak acid to determination the dissociation constant:

$$pH = pK_a + \log \frac{b}{(a - b)}$$

When concentration of both base and acid, at equivalent point  $a = b$  and at half equivalent point  $b = a/2$ ,  $(pH)_{1/2} = pK_a$  from  $(pH)_{1/2}$  calculate  $K_a$  value for acid.

**Procedure**

1. The electrode are washed in distilled water, dried with filter paper and washed again by successive immersions in the 0.1 M sodium hydroxide solution.
2. The electrodes are immersed in 25 ml of sodium hydroxide solution contained in a clean beaker.
3. This solution titrated with the 50 ml of 0.05 M acetic acid and the pH recorded after each addition of titrate, And stirring of the mixture before taking a pH reading.
4. On completion of the experiment the electrode is washed and immersed in distilled water.

**Calculation**

The flowing graphs are then plotted:-

- a) pH vs. volume of acetic aced added Fig.1 and b)  $\frac{\Delta (\text{pH})}{\Delta (v)}$  vs. volume of acetic acid added Fig.2

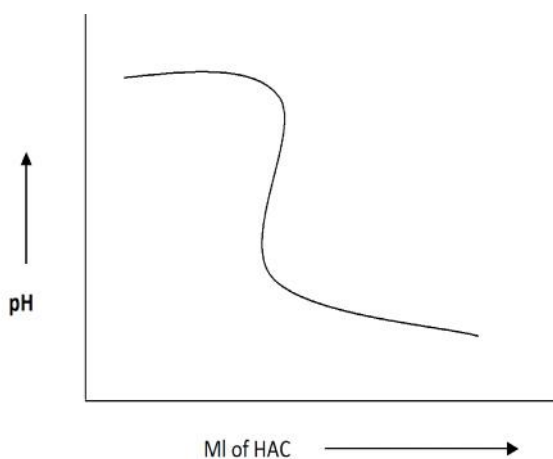


Fig.1

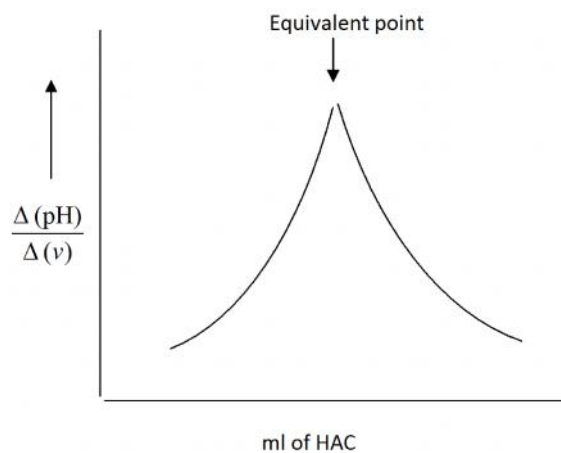


Fig.2

## Safety Guidelines in the Chemistry Laboratory

The chemistry laboratory is a safe place to work if you are aware of important safety rules and if you are careful. You must be responsible for your own safety and for the safety of others. The safety rules given here will protect you and others from harm in the lab. While carrying out procedures in any of the **ChemLabs**, notice the safety symbols and caution statements. The safety symbols are explained in the chart on the next page.












1. Always obtain your teacher's permission to begin a lab.
  2. Study the procedure. If you have questions, ask your teacher. Be sure you understand all safety symbols shown.
  3. Use the safety equipment provided for you. Goggles and a safety apron should be worn when any lab calls for using chemicals.
  4. When you are heating a test tube, always slant it so the mouth points away from you and others.
  5. Never eat or drink in the lab. Never inhale chemicals. Do not taste any substance or draw any material into your mouth.
  6. If you spill any chemical, wash it off immediately with water. Report the spill immediately to your teacher.
  7. Know the location and proper use of the fire extinguisher, safety shower, fire blanket, first aid kit, and fire alarm.
  8. Keep all materials away from open flames. Tie back long hair.
  9. If a fire should break out in the classroom, or if your clothing should catch fire, smother it with the fire blanket or a coat, or get under a safety shower. **NEVER RUN.**
  10. Report any accident or injury, no matter how small, to your teacher.
- Follow these procedures as you clean up your work area.
1. Turn off the water and gas. Disconnect electrical devices.
  2. Return materials to their places.
  3. Dispose of chemicals and other materials as directed by your teacher. Place broken glass and solid substances in the proper containers. Never discard materials in the sink.
  4. Clean your work area.
  5. Wash your hands thoroughly after working in the laboratory.






## First Aid in the Laboratory

Injury	Safe Response
Burns	Apply cold water. Call your teacher immediately.
Cuts and bruises	Stop any bleeding by applying direct pressure. Cover cuts with a clean dressing. Apply cold compresses to bruises. Call your teacher immediately.
Fainting	Leave the person lying down. Loosen any tight clothing and keep crowds away. Call your teacher immediately.
Foreign matter in eye	Flush with plenty of water. Use eyewash bottle or fountain.
Poisoning	Note the suspected poisoning agent and call your teacher immediately.
Any spills on skin	Flush with large amounts of water or use safety shower. Call your teacher immediately.

## Safety Symbols

These safety symbols are used in laboratory and investigations in this book to indicate possible hazards. Learn the meaning of each symbol and refer to this page often. *Remember to wash your hands thoroughly after completing lab procedures.*

SAFETY SYMBOLS	HAZARD	EXAMPLES	PRECAUTION	REMEDY
<b>DISPOSAL</b> 	Special disposal procedures need to be followed.	certain chemicals, living organisms	Do not dispose of these materials in the sink or trash can.	Dispose of wastes as directed by your teacher.
<b>BIOLOGICAL</b> 	Organisms or other biological materials that might be harmful to humans	bacteria, fungi, blood, unpreserved tissues, plant materials	Avoid skin contact with these materials. Wear mask or gloves.	Notify your teacher if you suspect contact with material. Wash hands thoroughly.
<b>EXTREME TEMPERATURE</b> 	Objects that can burn skin by being too cold or too hot	boiling liquids, hot plates, dry ice, liquid nitrogen	Use proper protection when handling.	Go to your teacher for first aid.
<b>SHARP OBJECT</b> 	Use of tools or glassware that can easily puncture or slice skin	razor blades, pins, scalpels, pointed tools, dissecting probes, broken glass	Practice common-sense behavior and follow guidelines for use of the tool.	Go to your teacher for first aid.
<b>FUME</b> 	Possible danger to respiratory tract from fumes	ammonia, acetone, nail polish remover, heated sulfur, moth balls	Make sure there is good ventilation. Never smell fumes directly. Wear a mask.	Leave foul area and notify your teacher immediately.
<b>ELECTRICAL</b> 	Possible danger from electrical shock or burn	improper grounding, liquid spills, short circuits, exposed wires	Double-check setup with teacher. Check condition of wires and apparatus.	Do not attempt to fix electrical problems. Notify your teacher immediately.
<b>IRRITANT</b> 	Substances that can irritate the skin or mucous membranes of the respiratory tract	pollen, moth balls, steel wool, fiberglass, potassium permanganate	Wear dust mask and gloves. Practice extra care when handling these materials.	Go to your teacher for first aid.
<b>CHEMICAL</b> 	Chemicals that can react with and destroy tissue and other materials	bleaches such as hydrogen peroxide; acids such as sulfuric acid, hydrochloric acid; bases such as ammonia, sodium hydroxide	Wear goggles, gloves, and an apron.	Immediately flush the affected area with water and notify your teacher.
<b>TOXIC</b> 	Substance may be poisonous if touched, inhaled, or swallowed.	mercury, many metal compounds, iodine, poinsettia plant parts	Follow your teacher's instructions.	Always wash hands thoroughly after use. Go to your teacher for first aid.
<b>FLAMMABLE</b> 	Open flame may ignite flammable chemicals, loose clothing, or hair.	alcohol, kerosene, potassium permanganate, hair, clothing	Avoid open flames and heat when using flammable chemicals.	Notify your teacher immediately. Use fire safety equipment if applicable.
<b>OPEN FLAME</b> 	Open flame in use, may cause fire.	hair, clothing, paper, synthetic materials	Tie back hair and loose clothing. Follow teacher's instructions on lighting and extinguishing flames.	Always wash hands thoroughly after use. Go to your teacher for first aid.

	<b>Eye Safety</b> Proper eye protection should be worn at all times by anyone performing or observing science activities.		<b>Clothing Protection</b> This symbol appears when substances could stain or burn clothing.		<b>Animal Safety</b> This symbol appears when safety of animals and students must be ensured.		<b>Radioactivity</b> This symbol appears when radioactive materials are used.		<b>Handwashing</b> After the lab, wash hands with soap and water before removing goggles
---	--	---	---	---	--	---	--	---	---



Physical Constants		
Quantity	Symbol	Value
Atomic mass unit	amu	$1.6605 \times 10^{-27}$
Avogadro's number	$N$	$6.022 \times 10^{23}$ particles/mole
Ideal gas constant	$R$	8.31 L·kPa/mol·K 0.0821 L·atm/mol·K 62.4 mm Hg·L/mol·K 62.4 torr·L/mol·K
Mass of an electron	$m_e$	$9.109 \times 10^{-31}$ kg $5.485799 \times 10^{-4}$ amu
Mass of a neutron	$m_n$	$1.67492 \times 10^{-27}$ kg 1.008665 amu
Mass of a proton	$m_p$	$1.6726 \times 10^{-27}$ kg 1.007276 amu
Molar volume of ideal gas at STP	$V$	22.414 L/mol
Normal boiling point of water	$T_b$	373.15 K 100.0°C
Normal freezing point of water	$T_f$	273.15 K 0.00°C
Planck's constant	$h$	$6.6260693 \times 10^{-34}$ J·s
Speed of light in a vacuum	$c$	$2.997925 \times 10^8$ m/s



## Alphabetical Table of the Elements

Element	Symbol	Atomic number	Atomic mass	Element	Symbol	Atomic number	Atomic mass
Actinium	Ac	89	227.027 8*	Mendelevium	Md	101	258.098 6*
Aluminum	Al	13	26.981 539	Mercury	Hg	80	200.59
Americium	Am	95	243.061 4*	Molybdenum	Mo	42	95.94
Antimony	Sb	51	121.757	Neodymium	Nd	60	144.24
Argon	Ar	18	39.948	Neon	Ne	10	20.179 7
Arsenic	As	33	74.921 59	Neptunium	Np	93	237.048 2
Astatine	At	85	209.987 1*	Nickel	Ni	28	58.6934
Barium	Ba	56	137.327	Niobium	Nb	41	92 .906 38
Berkelium	Bk	97	247.070 3*	Nitrogen	N	7	14. 006 74
Beryllium	Be	4	9.012 182	Nobelium	No	102	259.100 9*
Bismuth	Bi	83	208.980 37	Osmium	Os	76	190.2
Bohrium	Bh	107	262*	Oxygen	O	8	15.999 4
Boron	B	5	10.811	Palladium	Pd	46	106.42
Bromine	Br	35	79.904	Phosphorus	P	15	30.973 762
Cadmium	Cd	48	112.411	Platinum	Pt	78	195.08
Calcium	Ca	20	40.078	Plutonium	Pu	94	244.064 2*
Californium	Cf	98	251.079 6*	Polonium	Po	84	208.982 4*
Carbon	C	6	12.011	Potassium	K	19	39.098 3
Cerium	Ce	58	140.115	Praseodymium	Pr	59	140.907 65
Cesium	Cs	55	132.905 43	Promethium	Pm	61	144.912 8*
Chlorine	Cl	17	35.452 7	Protactinium	Pa	91	231.035 88
Chromium	Cr	24	51.996 1	Radium	Ra	88	226.025 4
Cobalt	Co	27	58.933 20	Radon	Rn	86	222.017 6*
Copper	Cu	29	63.546	Rhenium	Rc	75	186.207
Curium	Cm	96	247.070 3*	Rhodium	Rh	45	102.905 50
Dubnium	Db	105	262*	Rubidium	Rb	37	85.467 8
Dysprosium	Dy	66	162.50	Ruthenium	Ru	44	101.07
Einsteinium	Es	99	252.082 8*	Rutherfordium	Rf	104	261*
Erbium	Er	68	167.26	Samarium	Sm	62	150.36
Europium	Eu	63	151.965	Scandium	Sc	21	44.955 910
Fermium	Fm	100	257.095 1*	Seaborgium	Sg	106	263*
Fluorine	F	9	18.998 403 2	Selenium	Se	34	78.96
Francium	Fr	87	223.019 7*	Silicon	Si	14	28.085 5
Gadolinium	Gd	64	157.25	Silver	Ag	47	107.868 2
Gallium	Ga	31	69.723	Sodium	Na	11	22.989 768
Germanium	Ge	32	72.61	Strontium	Sr	38	87.62
Gold	Au	79	196.966 54	Sulfur	S	16	32.066
Hafnium	Hf	72	178.49	Tantalum	Ta	73	180.947 9
Hassium	Hs	108	265*	Technetium	Tc	43	97.907 2*
Helium	He	2	4.002 602	Tellurium	Te	52	127.60
Holmium	Ho	67	164.930 32	Terbium	Tb	65	158.925 34
Hydrogen	H	1	1.007 94	Thallium	Tl	81	204.383 3
Indium	In	49	114.82	Thorium	Th	90	232.038 1
Iodine	I	53	126.904 47	Thulium	Tm	69	168.934 21
Iridium	Ir	77	192.22	Tin	Sn	50	118.710
Iron	Fe	26	55.847	Titanium	Ti	22	47.88
Krypton	Kr	36	83.80	Tungsten	W	74	183.85
Lanthanum	La	57	138.905 5	Uranium	U	92	238.028 9
Lawrencium	Lr	103	260.105 4*	Vanadium	V	23	50.941 5
Lead	Pb	82	207.2	Xenon	Xe	54	131.29
Lithium	Li	3	6.941	Ytterbium	Yb	70	173.04
Lutetium	Lu	71	174.967	Yttrium	Y	39	88.905 85
Magnesium	Mg	12	24.305 0	Zinc	Zn	30	65.39
Manganese	Mn	25	54.938 05	Zirconium	Zr	40	91.224
Meitnerium	Mt	109	266*				

\* The mass of the isotope with the longest known half-life.

